
uVision Air

OS-1202 Instruction Sheet

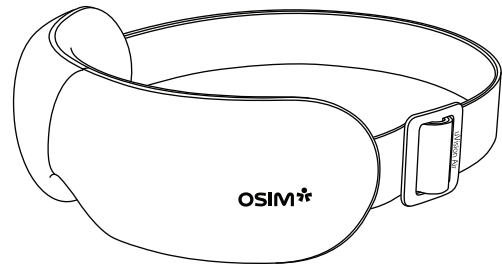Thank you for purchasing the OSIM uVision Air. Please read this instruction sheet thoroughly before using the uVision Air. It is your guide to a safe and efficient operation.

Safety

- This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities, or lack of experience or knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved.
- Children aged 7 and below shall not play with the product.
- The product must only be supplied at safety extra low voltage corresponding to the marking on appliance.
- The product has a heated surface. Persons insensitive to heat must be careful when using the product.
- Do not use the product if you are not feeling well.
- Do not store the product where it will be exposed to direct sunlight, excessive dust or humidity.
- Do not use the product in temperatures higher than 40°C.
- Do not wash the product or immerse it in water.
- Do not use the product in an environment with water.


- Keep the product away from pets and children.
- The product must be disconnected from the supply mains when removing the battery.
- This product contains batteries that are only replaceable by skilled persons.
- The battery must be removed from the product before it is scrapped.
- The battery is to be disposed of safely.
- This product contains batteries that are non-replaceable by user, do not attempt to disassemble or repair the product. Contact OSIM Service Center* for any examination, repair and battery removal.

*Please visit us at www.OSIM.com for information on the nearest service center.

⚠ Caution

- Only use the USB cable provided for charging.
- Do not leave the product unattended during charging.
- Do not use or rest on the product during charging.
- Product should be charged in an environment between 0 to 40°C.
- Please be mindful of your fingers when folding and unfolding the uVision Air.

Name of parts and operation

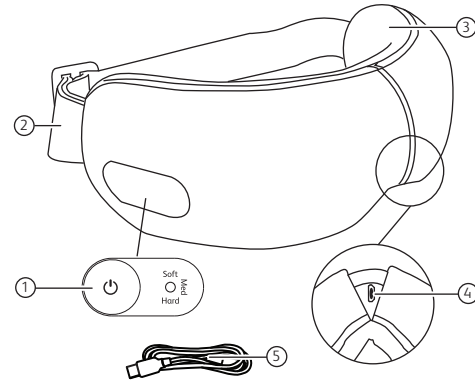
1. Power ON/OFF & Program selector switch
- Press and hold [⏻] for 2 seconds to switch product on/off. (Program will be at 'Med' when on)
- Press [⏻] to cycle through the 3 auto programs. (Soft, Med, Hard)
2. Adjustable, elastic head strap
3. Massage pad
4. Micro USB jack
5. USB cable

Charging the uVision Air

- Connect the USB cable to electronic devices such as laptop, USB wall power adaptor or portable charger.
- Insert the USB cable into the micro USB connector jack of the uVision Air.
- The product will not operate during charging mode.
- Each full charge (approximately 2-3 hours) lasts up to 240 min of usage (Usage time varies base on the program selected).

LED light indicator

- The LED indicator will turn ON and emit a blinking white light when battery level is running low.
- The LED indicator will turn ON and emit a blinking white light while charging.
- The LED indicator will turn ON and emit a steady white light when the battery is fully charged.

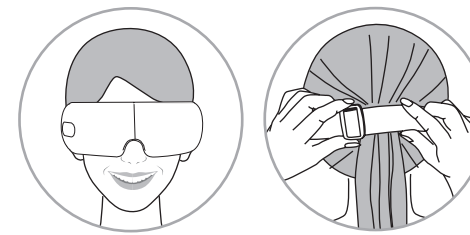Fitting on the uVision Air

- The uVision Air is to be used as shown above only. Covering your eyes completely.
- Tighten or loosening the product by adjusting the strap until you feel comfortable. To experience a better massage, it is recommended that the product is strapped securely on your face.

Product care and storage

- Do not stack heavy items on the product.
- Before cleaning, ensure that the product is turned off and disconnected from the power source.
- Use a soft, dry lint-free cloth for daily cleaning. Do not use volatile liquids such as benzene, thinner or gasoline for cleaning.

Specifications

Model	OSIM uVision Air
Model Number	OS-1202
USB Input	5 V 1A
Dimension	21 (L) X 7.5 (W) X 4 (H) cm
Net Weight	Approx. 300g
Gross Weight	Approx. 770g

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